
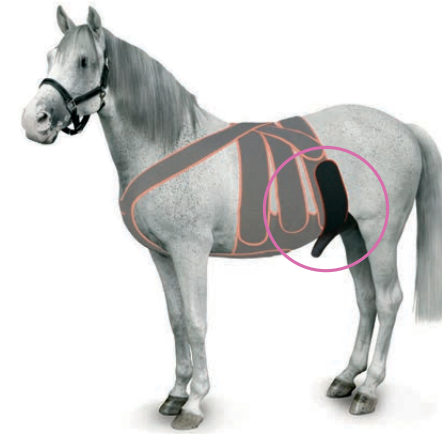


# KRUUSE Equine Recovery & Rehabilitation Bandages

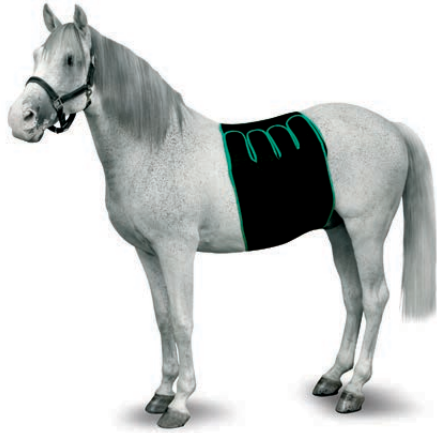

## KRUUSE Recovery Bandaging

**Indication:** For immediate post-operative protection for the first hours after the surgery.

A 3-dimensional elastic wrap with wide foam endings covered in velour to form a strong adhesive to TRIHOOK.

Interior pocket lined with a monofilm pad 1cm thick, acts as an air cushion for covering the lower abdomen. The KRUUSE Recovery Bandage allows for the horse to stand up easier and provides full protection in the case of a violent recovery.

Provet Code	Description	
EQUI K R1	Foal Recovery Bandaging S: 70 - 89cm	●
EQUI K R2	Foal Recovery Bandaging M: 90 - 109cm	●
EQUI K R3	Foal Recovery Bandaging L: 110 - 130cm	●
EQUI K R4	Equine Recovery Bandaging S: 130 - 170cm	●
EQUI K R5	Equine Recovery Bandaging M: 170 - 200cm	●
EQUI K R6	Equine Recovery Bandaging L: 200 - 240cm	●

## KRUUSE Post Colic Surgery Kit

**Indication:** Abdominal surgery.

An abdominal wrap with a removable and washable interior pocket. The mesh inside the interior pocket may be slit to create an opening for inserting and replacing absorbent drainage pads.

Fit in less than three minutes the perfect post-operative bandaging.

External: Foam covered in velour that is adhesive to TRIHOOK (an extra strong adhesive strip).

Internal: 100% skin-friendly cotton.

Provet Code	Description	
POST COL1	Foal Post Colic Surgery Kit S: 70 - 89cm	●
POST COL2	Foal Post Colic Surgery Kit M: 90 - 109cm	●
POST COL3	Foal Post Colic Surgery Kit L: 110 - 130cm	●
POST COL4	Equine Post Colic Surgery Kit S: 130 - 170cm	●
POST COL5	Equine Post Colic Surgery Kit M: 170 - 200cm	●
EQUI P 5	Equine Post Colic Surgery Kit L: 200 - 240cm	●
POST COL 7	Colic Recovery Replacement Pockets	

## KRUUSE Equine Scrotum Strap

**This strap is an addition to the Post Colic Surgery Kit and only for male horses.**

The strap has three adjustment points. The goal is to adjust the small hollow just in front of the penis to avoid the horse urinating inside the abdominal wrap.

**Use:** Fit the KRUUSE Scrotum Strap in the middle just in front of the penis and fasten on the dorsal strap.

Provet Code	Description	
EQUI K S1	Equine Scrotum Strap S	●
EQUI K S2	Equine Scrotum Strap M	●
EQUI K S3	Equine Scrotum Strap L	●

**Note:** Choose same size as post colic surgery kit

Scan QR codes above to view product videos

## Provet Branch Contacts

P 0800 776 838 F 0800 776 839



A covetrus Company

### Auckland

CServiceAuckland@provet.co.nz


### Christchurch

CServiceChristchurch@provet.co.nz

### Palmerston North

CServicePalmerstonN@provet.co.nz

## How to measure the size

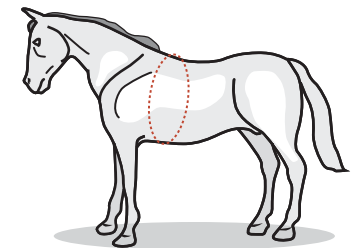

Size Colour Code:

● Small ● Medium ● Large